

# DIGITAL BUTYRO REFRACTOMETER


# PR-BUTYRO

Cat.No.3454


New ATAGO Digital Model PR-BUTYRO perfect for Butter, Ghee & Edible Oil industry for checking purity / adulteration.

## Features

- Calibration with Standard liquid
- Large, easy to read display
- External Light Interference (ELI)\*
- Display of both Butyro & RI scale
- Temperature Compensated Scale at 40°C
- Automatic & Accurate (Accuracy : Butyro±0.5)

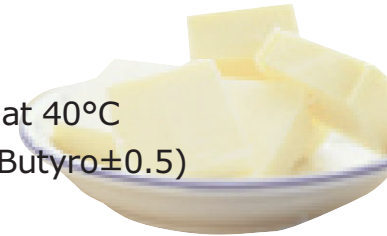

Dual scale of Butyro & Refractive Index !

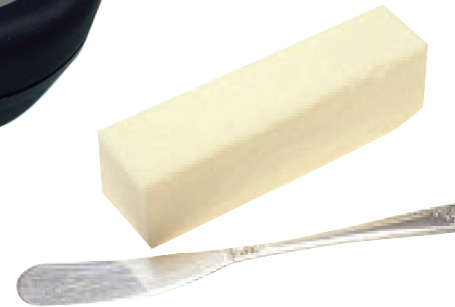

## Specifications

Measurement Range (Automatic Temperature Compensation)	Butyro 30.0 to 90.0 RI 1.4450 to 1.4850 (Converted at 40°C)
Minimum Indication	Butyro 0.1 RI 0.0001
Measurement Accuracy	Butyro ±0.5 (at 40°C) RI ±0.0003 (at 40°C)
Measurement Temperature	10 to 50°C (Automatic Temperature Compensation)
Ambient Temperature	10 to 40°C
Sample Volume	0.1ml or more
Measurement Time	3 seconds (approximately)
Power Supply	006P dry battery (9V)
International Protection Class	IP 64 Water Resistant
Dimensions and Weight	17(W)×9(D)×4(H)cm, 300g (Main unit only)

### \*External Light Interference (ELI)

When intense light penetrates the prism of a digital refractometer, the light waves may interfere with the sensor resulting in inaccurate readings.

To ensure accurate measurement results each time, the PR-BUTYRO is programmed with the ELI function to display the warning message "nnn" when intense direct light is detected. Forming a habit of shading the sample stage with your hand and re-pressing the START/OFF Switch (when the warning message from the ELI function is displayed) will ensure accurate measurement results each time.

All ATAGO refractometers are designed and manufactured in Japan.



Headquarters: The Front Tower Shiba Koen,  
23rd Floor 2-6-3 Shiba-koen, Minato-ku,  
Tokyo 105-0011, Japan

TEL : 81-3-3431-1943 FAX : 81-3-3431-1945



ATAGO products comply with HACCP,GMP, and GLP system standards.

<http://www.atago.net/> [overseas@atago.net](mailto:overseas@atago.net)

**ATAGO U.S.A., Inc.**

TEL : 1-425-637-2107

customerservice@atago-usa.com

**ATAGO INDIA Instruments Pvt. Ltd.**

TEL : 91-22-2833-8038 / 8076

customerservice@atago-india.com

**ATAGO (THAILAND) Co., Ltd.**

TEL : 662-982-8718-9

customerservice@atago-thailand.com

**ATAGO BRASIL Ltda.**

TEL : 55 16 3916-6000

customerservice@atago-brasil.com

**ATAGO ITALIA s.r.l.**

TEL : 39 2 3657267

customerservice@atago-italia.com

**ATAGO CHINA Guangzhou Co., Ltd.**

TEL : 86 20-38106065

info@atago-china.com

\* Specifications and appearance are subject to change without notice.